

CASE STUDY

Purity and identity characterization of AAV capsid particles by LC-MS methods

Wendi A. Hale¹, Dominique Garceau², Tristan Cano², Caitlin Jaeger², Roy Hegedus², William Hermans², Norman Garceau², Christopher M. Colangelo¹

¹Agilent, Lexington, MA; ²Curia, Albany, NY

Introduction

Adeno-associated viruses (AAVs) are the main viral vectors for gene therapy and have been successful in treating inherited retinal diseases and spinal muscular atrophy. AAV is composed of an icosahedral protein shell with a single-stranded genome of approximately 4.7 kb. The intact AAVs act as a vehicle to protect and deliver oligonucleotide therapeutics. There are 11 known serotypes that transduce different cell types, allowing for increased selectivity for therapeutics. As AAVs continue to be explored as therapeutic delivery platforms, it is vital to ensure that all the critical quality attributes (CQAs) of the therapeutic product are maintained.

Characterizing viral capsid proteins yields several challenges. The protein shell is composed of three capsid proteins, VP1, VP2 and VP3, that assemble into a 3.9 megadalton structure in a ratio of 1:1:10 with 60 capsids per virion. In addition to the low molar ratios of VP1 and VP2, all three proteins have overlapping sequences at the C-terminus. Traditionally, SDS-PAGE is used to establish the molecular weight of the

capsid proteins, however, this technique provides an approximate molecular weight and may not be able to distinguish between different serotypes. Mass spectrometry is a promising method to overcome these challenges and determine CQAs of the capsid proteins.

Experimental

MATERIALS:

AAV8 was produced by Lake Pharma (Worcester, MA). Molecular weight cutoff filters and (tris(2-carboxyethyl)phosphine) (TCEP) were purchased from Millipore Sigma. Trypsin and rAsp-N were purchased from Promega.

SAMPLE PREPARATION:

For intact analysis, AAVs underwent a buffer exchange three times at 10,000 g with a 10 kDa molecular weight filter. The buffer contained 5 mM TC EP, 20% H₂O and 80% acetonitrile with 0.1% formic acid (v/v). After the sample was collected, it was incubated at room temperature prior to injection. For peptide mapping, the AAVs underwent denaturation, reduction, alkylation and digestion. Enzymes utilized in this experiment were trypsin and rAsp-N.

Experimental

LC/MS ANALYSIS:

LC/MS analysis was performed on a 1290 Infinity II LC coupled to a 6545XT AdvanceBio LC/Q-TOF system with a dual Agilent Jet Stream source. Agilent MassHunter® Acquisition (B.09.00) workstation software with the large molecule SWARM autotune feature was used for intact analysis. The instrument was further calibrated and operated in standard mass mode. The iterative MS/MS feature was used for the peptide mapping workflow. All MS data was processed with Agilent MassHunter BioConfirm 10.0 software.

Figure 1: 6545XT AdvanceBio LC/Q-TOF



Table 1: Column and LC Conditions for Both Analyses

| | Intact Analysis | Peptide Mapping Analysis |
|--------------------|---|---|
| Column | Zorbax® Dipheryl RRHD 300Å, 2.1 x150 mm, 1.8 µm | AdvanceBio® Peptide Mapping, 2.1 x 150 mm, 2.7 µm |
| Flow Rate | 0.4 ml/min | 0.4 ml/min |
| Injection Volume | 20 µL | 40 µL |
| Column Temperature | 60°C | 60°C |

Results and discussion

INTACT ANALYSIS OF AAV8 ON THE 6545XT ADVANCEBIO LC/Q-TOF

The spectral clarity provided by the improved vacuum on the 6545XT AdvanceBio LC/Q-TOF in combination with the large molecule SWARM autotune feature show all three viral capsid proteins with their post-translational modifications (PTMs) with high mass accuracy, under 10 ppm for all proteoforms. While it is not shown, VP1 and VP2 are chromatographically separated which can be challenging. While mass spectrometry can separate these proteins by mass, having chromatographic separation allows for less ion suppression of these two low abundant proteins.

Figure 2: VP1 raw and deconvoluted mass spectra. Three phosphorylation sites were detected on VP1 with less than 10 ppm error. The accurate mass data confirmed that VP1 is missing its N-terminal amino acid residue and that the new N-terminus is acetylated.

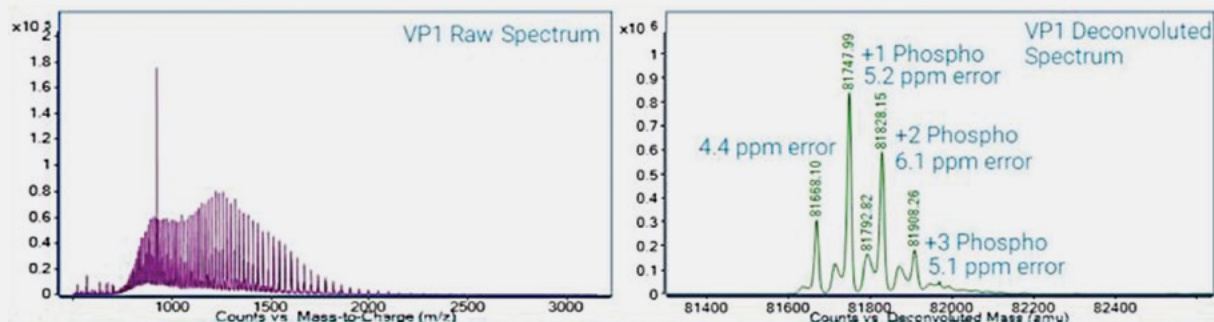


Figure 3: VP2 raw and deconvoluted mass spectra. The accurate mass data confirms at least two phosphorylation sites on VP2.


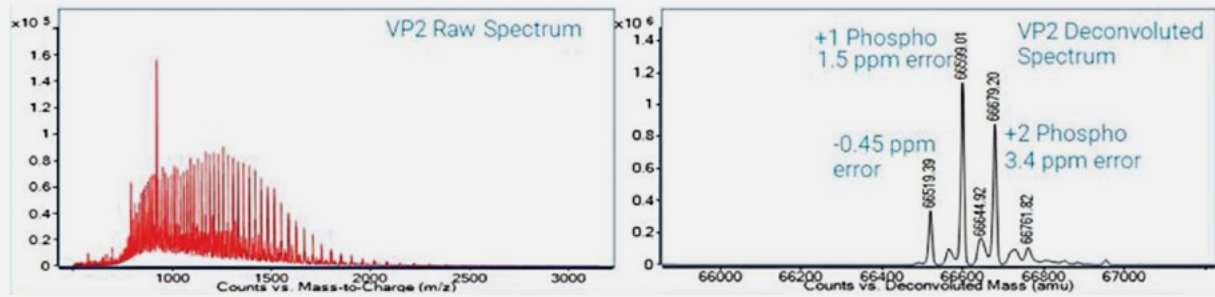


Figure 4: VP3 raw and deconvoluted mass spectra. The unmodified form of VP3 is mostly chromatographically separated from acetylated VP3. While VP1 was fully acetylated, about 70% of VP3 was acetylated.

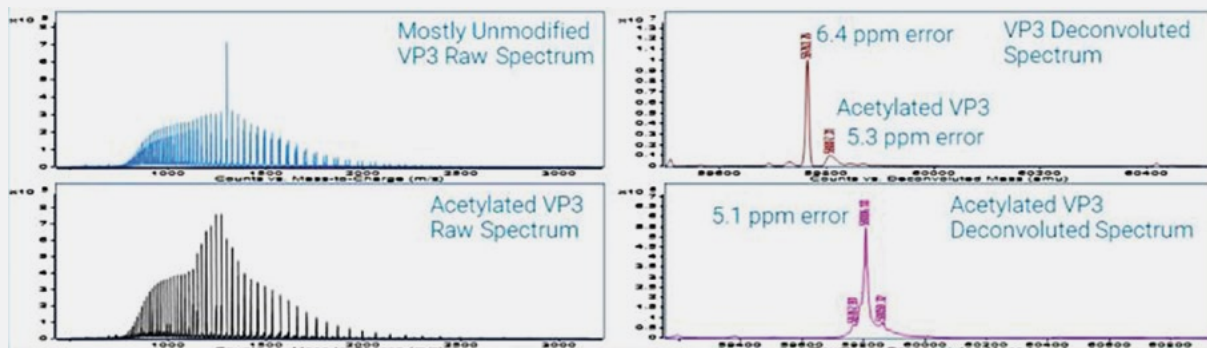

Results and discussion

PEPTIDE MAPPING OF VP1, VP2, AND VP3 ON THE 6545XT ADVANCEBIO LC/Q-TOF:

Peptide mapping of biotherapeutics is an essential method to determine protein sequence and post-translational modifications, required by the ICH, FDA and other regulatory agencies. As of January 2020, the FDA recommends providing information regarding primary and secondary structure including PTMs for human gene therapy drug substances. Peptide mapping with the iterative MS/MS feature excludes peptides from all previous runs for isolation and fragmentation, allowing for selection and detection of low abundant peptides. In addition, BioConfirm 10.0 allows for multiple runs to be selected to provide a total sequence coverage. This feature is useful for combining results from iterative MS/MS runs as well as using multiple enzymes.

Figure 5. Example of MS/MS spectrum confirming site-specific phosphorylation of a serine residue.

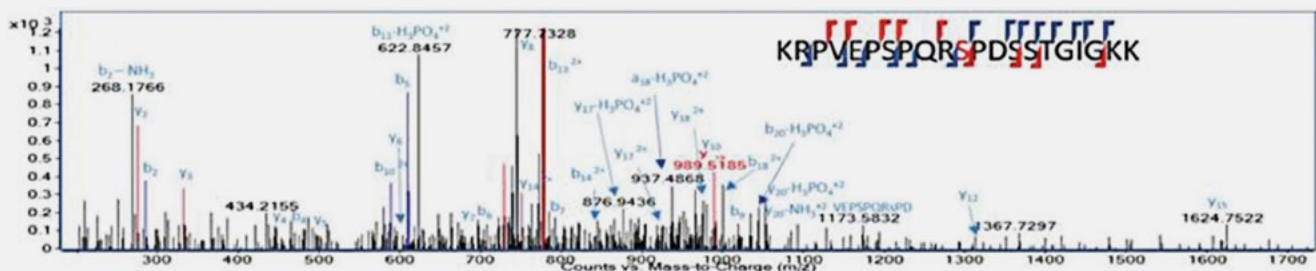


Figure 6. Peptide mapping results of each of the viral capsid proteins. Each protein has 100% or nearly 100% sequence coverage with MS/MS confirmation.


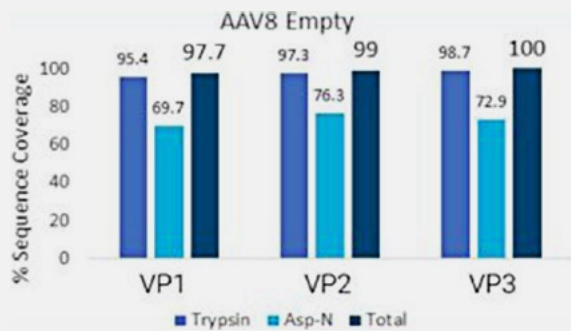


Figure 7. Examples of common PTMs. Low levels of deamidation and oxidation indicate this protein has not degraded.

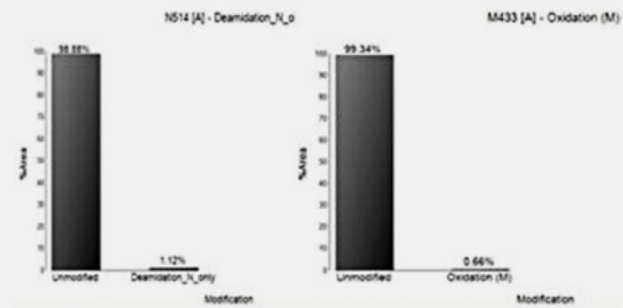

Conclusions

- All three capsid proteins were chromatographically separated.
- The intact protein data has clear spectra and all proteoforms have less than 10 ppm error.
- The peptide mapping data gives between 97.7-100% sequence coverage with MS/MS confirmation.
- Site specific phosphorylation was localized.
- Relative quantitation of PTMs was performed.

References

- 1 Dalkara, D. et al. *Sci. Transl. Med* 2013. 5(189), 189ra76-189ra76.
- 2 Moore, N. A. et al. *Expert Opin. Biol. Ther.* 2018. 18(1), 37-49.
- 3 Wu, Z. et al. *J. Mol. Ther.* 2006. 14(3), 316-327.
- 4 Jin, X. et al. *Hum. Gene Ther. Methods* 2017. 28(5), 255-267.
- 5 Giles, AR. et al. *Mol. Ther.* 2018. 26(12), 2848-2862.

For Research Use Only. Not for use in diagnostic procedures.

ABOUT CURIA

Curia is a Contract Development and Manufacturing Organization with over 30 years of experience, an integrated network of 29 global sites and over 3,500 employees partnering with customers to make treatments broadly accessible to patients. Our biologics and small molecule offering spans discovery through commercialization, with integrated regulatory and analytical capabilities. Our scientific and process experts and state-of-the-art facilities deliver best-in-class experience across drug substance and drug product manufacturing. From curiosity to cure, we deliver every step to accelerate and sustain life-changing therapeutics. ***Learn more at curiaglobal.com***

Solutions developed by Curia

CONTACT US

www.curiaglobal.com

